

OPERATION MANUAL

MEGAFLOW Peristaltic Pump – Tubing Pump



LAMBDA MEGAFLOW Peristaltic Pump – Tubing Pump

The LAMBDA peristaltic pumps have been developed for continuous cultures as the result of over twenty years of laboratory experience and involved the systematic elimination of the imperfections found in other pumps on the market.

The successful design and well-proved mechanics of the LAMBDA PRECIFLOW pump has been extended by flow rate programming. Up to 99 steps of time and flow rate can be easily programmed, thus allowing the creation of any desired flow rate profile. The maximum flow rate has been increased to up to 60,000 ml/hour. Until now, it was not possible to produce peristaltic pumps with such a high flow rate in such a small instrument casing.

- Most compact peristaltic pump of this type on the market
- Flow rates from 0.2 to 60,000 ml/hour
- Large digital speed setting range from 0 to 999
- Greatly extended tubing life and decreased pulsation
- Flow rate programming (up to 99 steps) and automatic switch-on and -off without using any timer
- Extensive remote controls
- Very economic in use
- Access to reaction kinetics by using the pump-flow INTEGRATOR
- Low voltage plug-in power supply for maximum safety
- RS-485 interface (optional)
- Control software PNet (optional)

LAMBDA Laboratory Instruments

is developer and producer of special laboratory instruments mainly for biotechnology, microbiology, food and agricultural, chemical and pharmaceutical research and development as well as for general laboratory and research applications.

LAMBDA MINIFOR – highly innovative and compact fermenter and bioreactor system for laboratory scale fermentation and cell cultures

LAMBDA OMNICOLL – fraction collector-sampler for unlimited number of fractions

LAMBDA PRECIFLOW, MULTIFLOW, HIFLOW and MAXIFLOW peristaltic pumps – reliable, precise and extremely compact

LAMBDA SAFETY POWDER DOSER – allows automatic feeding of powders without spoon. Safe operation with hazardous material (GLP)

LAMBDA VIT-FIT polyvalent syringe pump with extremely robust mechanics – programmable infusion and filling from micro syringes to large volume syringes of 150 ml without adapter

LAMBDA MASSFLOW – precise gas flow measurement and control with data acquisition option

LAMBDA PUMP-FLOW INTEGRATOR – with LAMBDA pumps and doser allows the visualization and recording of the pumped volume

Table of contents

1	Setting up the Peristaltic Pump	3
1.1	Tubing Insertion	3
1.2	ON/OFF button	3
1.3	Setting up the flow rate	3
1.4	Choosing flow direction	4
1.5	Fast filling or emptying the line	5
1.6	FAS / SLO mode	5
2	Programming the peristaltic pump	6
3	Peristaltic pump flow calibration	10
3.1	Volumetric calibration of the peristaltic pump flow	10
3.2	Peristaltic pump flow calibration by weight	11
4	Remote controls	12
4.1	ON/OFF remote control	12
4.2	Remote control of the pump speed	12
4.3	PC control	13
5	Recommendations	14
6	For your safety	14
7	Construction advantages of the pump	15
8	Utilization of LAMBDA Peristaltic Pumps	15
9	Technical specifications	16
9.1	General specification	16
9.2	Remote control (Inputs/outputs)	17
9.3	Input (12 V DC)	17
10	Guarantee	17
11	ACCESSORIES AND SPARE PARTS	18
11.1	Pump flow integrator (Art. No. 4803)	18
11.2	Pnet control software for peristaltic and syringe pumps, DOSER or MASSFLOW (Art. No. 6600)	18
11.3	List of accessories and spare parts	19
12	Appendix	20
12.1	RS communication protocol for LAMBDA PERISTALIC and SYRINGE PUMPS, powder DOSER and MASSFLOW gas flow regulator	20
12.2	Examples	21
12.3	How to set the MEGAFLOW peristaltic pump address?	21
12.4	RS-connection scheme	21
12.5	RS communication protocol for the on-board INTEGRATOR (optional)	22

1 SETTING UP THE PERISTALTIC PUMP

1.1 Tubing Insertion

A short video of the peristaltic pump installation can be viewed online: <http://www.lambda-instruments.com/?pages=video#peristaltic>

- Plug the connector of the power supply into the corresponding socket on the rear side of the peristaltic pump and secure it in place.
- Plug the power supply into the AC mains (230V/50-60 Hz).
- Remove the clear PVC cover by rotating it, to the left or to the right.
- Press the **ON/OFF** button and select the sense of rotation using the button ◀|▶, so that the rotor will turn in a clockwise direction.
- Set the pump speed to about **100** (sense of rotation ▶) by pressing the corresponding buttons ▲▲▲ below the LED display.
- Press the tubing into the back slot on the top of the peristaltic pump. Thin tubing should be pushed completely to the bottom of the slot. Guide the tubing around the outside of the slowly turning plastic bearing towards the front slot. Then, press the tubing into the front slot to secure it.
- Place the transparent cover on the top of the pump and turn it, so that the steel ball embedded in the left-front corner fits into the corresponding notch on the PVC cover. When the cover is gently pressed, the bearings will automatically adjust themselves.



Care should be taken not to get the fingers clamped by the rollers, while inserting the tubing on the pump head

1.2 ON/OFF button

By pressing the **ON/OFF** button, the peristaltic pump is switched on or off. The internal memory will show the last used speed and flow direction setting.



Figure 1-1 Pressing **ON/OFF** button to switch **ON** the peristaltic pump



Figure 1-2 The peristaltic pump will run at the set speed and flow direction

1.3 Setting up the flow rate

The flow rates delivered by peristaltic pumps depend on the internal diameter of the tubing and the pump rotation speed. The speed of rotation can be selected by the control buttons ▲▲▲ under the LED display.

The MEGAFLOW tubing pump has been constructed for tubing with an internal diameter ranging from 1 to 8 mm with a tubing wall thickness of approximately 2 mm. The best results have been obtained with silicone tubing, but tubing made from other materials with similar elasticity can also be used.

The flow rates of the LAMBDA MEGAFLOW peristaltic pump can be varied over a large range from 0.2 to 60,000 ml/hour.



Figure 1-3 The ▲▲▲ buttons can be used to set the desired flow rate.



Figure 1-4 Each ▲▲▲ arrow can be used to change the corresponding digit on the flow rate.

With the control buttons ▲▲▲ below the LED display the motor speed is selected. The speed setting from 0 to 999 corresponds to the velocity of the movement of the motor. The best way to correlate the flow rate obtained with the respective tubing is to make a preliminary calibration, in which the pump is allowed to pump the liquid over a certain time with a selected speed setting (e.g. for 1 minute with speed setting 500). Then, the volume (*refer section 3.1*) or weight (*refer section 3.2*) of the pumped sample are measured. Using this information, the speed setting corresponding to the desired flow rate can be calculated easily (rule of three).

1.4 Choosing flow direction

The direction of the pump rotation can be selected by the ◀|▶ button, clockwise or anti-clockwise. The corresponding direction LED diode will be ON.



If possible use the clockwise rotation of the tubing pump. This results in lower friction and the pressure of the liquid (approx. 0.18 MPa). If a higher pressure is required (up to 0.2 MPa), use the counter-clockwise rotation.



Figure 1-5 LED illumination in the direction ◀| represents the clockwise rotation.



Figure 1-6 LED illumination in the direction |▶ represents the anticlockwise rotation.

1.5 Fast filling or emptying the line

If the direction ◀|▶ button is pressed continuously for about 2 seconds, the pump will rotate at a maximum speed in the direction of rotation indicated by the LED.

After releasing this direction button, the pump will stop pumping. This factor can be used for filling tubing before starting or for emptying the tubing line at the end of the operation.

This “**HOLD=MAX**” function can be used even if the **ON/OFF** button has not been pressed.

1.6 FAS / SLO mode

The MEGAFLOW peristaltic pumps allow to pump the liquid at two speed settings: **SLO.** = slow and **FAS** = fast.

In the standard mode (**FAS**), the flow rate of the speed setting 0 to 999 ranges from 0 to 60,000 ml/hour.

In the slow mode (**SLO.**), the flow rate of the speed setting 0 to 999 ranges from 0 to 12,000 ml/hour. The slow mode is indicated by the dot in the last position on the digital display.

The fast filling function ◀|▶ (**HOLD = MAX**) is the same for both speed ranges.



Figure 1-7 Plug-in the power supply and at the same time press the **RUN** button on the control panel.



Figure 1-8 The display will show the default mode **FAS** or the last used mode (**FAS / SLO**).



Figure 1-9 Press the **RUN** button, to select the desired speed mode (**FAS / SLO**).



Figure 1-10 Confirm and save the chosen speed mode by pressing the **ON/OFF** button on the control panel.

Remark: During the RS-communication, only the value shown on the display is transferred, but not the speed range (mode). If the PUMP-FLOW INTEGRATOR is activated (optional), a single integrator step is independent of the set speed range (mode).

2 PROGRAMMING THE PERISTALTIC PUMP

A short video about programming the LAMBDA Peristaltic Pumps can be found at: <http://lambda-instruments.com/?pages=video#peristaltic>.

Up to 99 pairs of time and speed settings (flow rates) may be programmed in a simple way.

The programming mode is accessed by simultaneously pressing the buttons **REMOTE** and **RUN** until the indication “**PGM**” appears on the display and both directions LEDs (◀▶) are illuminated.



Figure 2-1 Press the **REMOTE** and **RUN** simultaneously until the indication “**PGM**” appears on the display.



Figure 2-2 Continuous pressing of **REMOTE** and **RUN** button even after the indication of “**PGM**” ends up with the “**cLE**” indication.

Remark: If you repeat this simultaneous pressing of the **REMOTE** and **RUN** buttons, the memory will be cleared and the indication “**cLE**” will appear on the display. To enter the programming mode again, press the **REMOTE** and **RUN** buttons again until “**PGM**” appears.



Figure 2-3 Press the **ON/OFF** button. The indication “**F01**” will appear for a short time on the display indicating that you can select the first flow rate (speed setting) value.



Figure 2-4 Set the desired flow rate value for the first program step by pressing the buttons ▲▲▲ below the display (from 0 to 999, corresponding to 0 to 100% of the motor rotation speed).



Figure 2-5 Select the direction of the pump rotation using the direction button ◀▶.



Figure 2-6 Press the **ON/OFF** button. The indication “t01” will appear for a few seconds on the display indicating that the time period of the first step in minutes, can be programmed.



Figure 2-7 Select the desired time period for the first program step by pressing the buttons ▲▲▲ below the display (from 0 to 999 minutes or 00.0 to 99.9 minutes).



Figure 2-8 By pressing the ◀▶ button, the time resolution can be set in minutes or 0.1 minutes. In the 0.1 minute time resolution a dot is displayed, e.g. “00.2”. The time resolution can be set individually for each program step.



Figure 2-9 Press the **ON/OFF** button. The indication “F02” will briefly appear on the display.



Figure 2-10 Enter the desired flow rate for the second program step.



Figure 2-11 If the ◀▶ button is not used to change the direction, the default pump rotation for the set speed is clockwise, ◀ LED ON.



Figure 2-12 Press the ON/OFF button again, the symbol “t02” will briefly appear on the display.



Figure 2-13 Set the time of the second program step.



Figure 2-14 Press ◀▶ button, to set the time resolution.

In a similar way up to 99 program steps can be entered.



Figure 2-15 After having entered the time of the last step, press the ON/OFF button.



Figure 2-16 The flow rate (000) of the next step which will not be programmed appears on the display.

Remark: It is not possible to end the program after programming the time data.

The direction LEDs indicate if you are programming speed or time:

- One direction LED is on: Programming of the flow rate (in the direction indicated by the LED)
- Both directions LED are off: Programming of the time



Figure 2-17 Press the **Remote** and **Run** button simultaneously and you will see the indication “c01” on the display. This indicates that the program will be executed only once and the peristaltic pump will stop afterwards.



Figure 2-18 If you wish to repeat the same program 3 times, increase the cycle number to “c03” by pressing the buttons $\wedge \wedge \wedge$ below the display (from 0 to 99 cycles). The program can be repeated up to 99 times, indicated by “c99”. If 0 is entered for the cycle number “c00”, the program will run continuously (infinite loop).



Figure 2-19 Press the **ON/OFF** button until the indication “End” appears on the display to confirm and save the program.



Figure 2-20 To start the program, press the **RUN** button. The **RUN** and **ON/OFF** LEDs are on. To stop the running program definitively, press the **RUN** button. The **RUN** and **ON/OFF** LEDs are off

It is possible to stop the pump by pressing the ON/OFF button, to change the direction and the rotation speed during any running program step. This facilitates urgent replacement of tubing or allows reaction in other emergency situations.

Remark: Do not forget to restore the right direction of the rotation and to switch the pump on again (by pressing the ON/OFF button) after you have finished your intervention.

The **time basis** in the microprocessor is **not stopped** during this intervention, so that the total time of the running steps and of the whole program will not be affected. **When the program step time has elapsed, the pump will automatically go on with the next program step.** Thus, the program is not modified by this emergency intervention.

It is possible to review the program by proceeding in the same way as during programming but without modifying it.

3 PERISTALTIC PUMP FLOW CALIBRATION



If you need to know the amount of liquid added without a balance or if you want to measure the flow rate, then you will need the calibration of your LAMBDA pump speed / flow rate.



While calibrating the pump flow, make sure that the liquid completely fills the tubing and reaches its free end. Only then initiate the measuring process, for precise calibration.

The calibration of the pump flow rate with speed can be done to know the amount of the liquid pumped.

A short video of peristaltic pump flow calibration can be found at <http://lambda-instruments.com/?pages=video#peristaltic>

The MEGAFLOW peristaltic pumps can be programmed for the calibration, for example: speed 300 and time 1 minute.

3.1 Volumetric calibration of the peristaltic pump flow

In the volumetric calibration, the amount of liquid pumped at a particular speed for a minute is calculated. Program the pump with flow rate say 600 and time of about 1 minute.

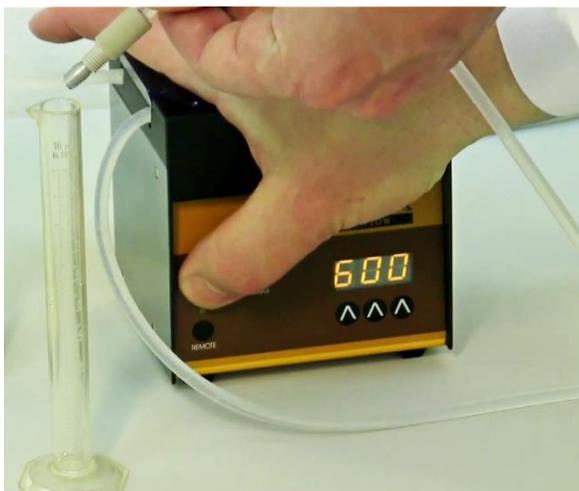


Figure 3-1 Turn on the pump and have the other end of the tubing ready near the measuring cylinder. Carefully collect the liquid being pumped in the measuring cylinder for about 60 seconds.

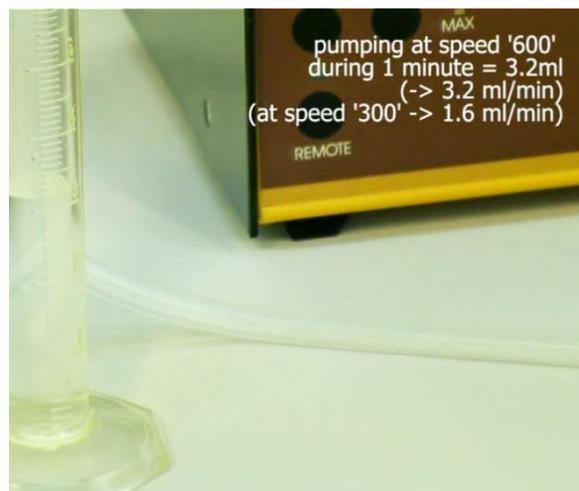


Figure 3-2 Measure the volume of liquid collected in a minute.

For example: At the pumping speed of 300 with the MEGAFLOW pump, the volume of the liquid collected is 10 ml/min. With this calculation, the volume of liquid collected at other pump speed can be found.

3.2 Peristaltic pump flow calibration by weight

In the calibration by weight, the weight of the liquid pump in a minute was calculated. Program the MEGAFLOW peristaltic pump for the calibration by weight.



Figure 3-3 Measure the weight of the empty beaker using a sensitive weighing scale. For the exact measurement of the liquid collected, tare the weighing scale with the measuring beaker (Example: 0.000 g).

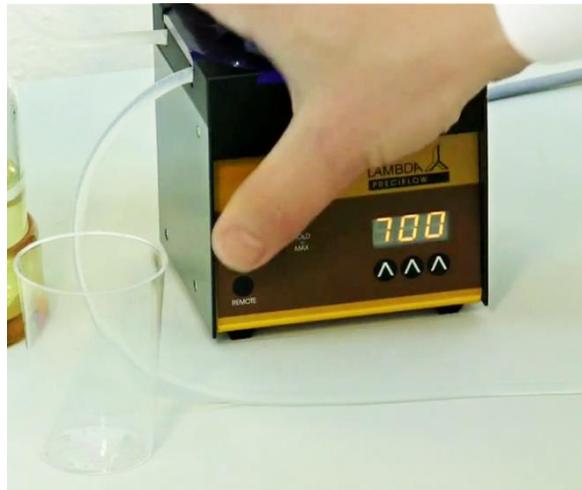


Figure 3-4 Connect the one end of the tubing to the liquid source and press ON/OFF or RUN button of the pump to start the program.

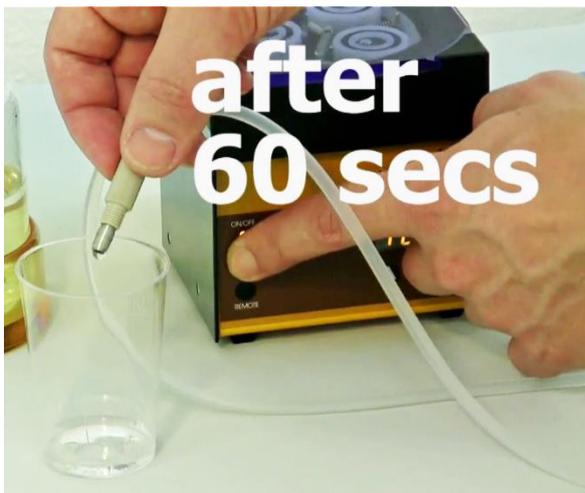


Figure 3-5 Carefully collect the liquid being pumped into the measuring beaker for about 60 seconds. Pump gets switched off exactly after 60 seconds as programmed.

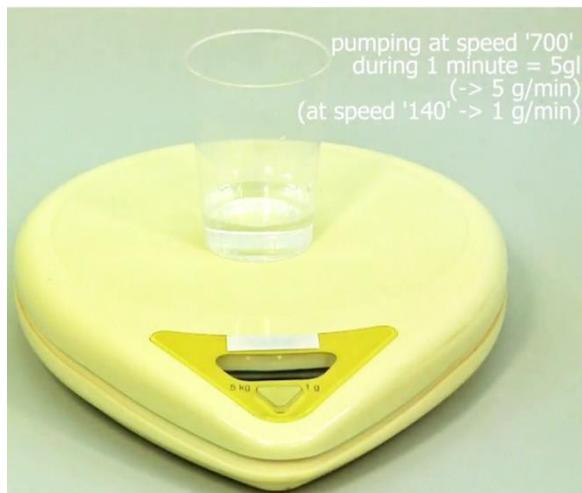


Figure 3-6 Weigh the measuring beaker with the liquid been collected in the past 60 seconds.

For example: At the pumping speed of 200 with the MEGAFLOW pump, the weight of the liquid collected is 8 g/min. With this calculation, the weight of liquid collected at other pump speed can be found.

4 REMOTE CONTROLS

4.1 ON/OFF remote control

By interlinking the contacts no. 4 and 5 of the socket at the rear of the pump (see *figure 9–1 and section 9.2*), the pump will be blocked and the ON/OFF LED will be switched off.

The same effect will be obtained by applying a voltage of 3 to 12 V DC to the contact no. 5 (0 line must be connected to contact no. 3).

Remark: In some cases a reversed logic for the remote control might be desired. Please contact us in this case.

4.2 Remote control of the pump speed

The LAMBDA Peristaltic pumps can be controlled over the whole speed range by an external signal (0 - 10 V DC, option 0-20 or 4-20 mA). The plus pole of the signal is connected to the contact no.1, 0 line to the contact no.3.



Figure 4-1 For remote control, connect the 8-poles Pump remote control cable into the 'Remote' socket at the rear of the Peristaltic Pump.



Figure 4-2 Press the button REMOTE on the front panel. The corresponding LED diode will go on and the display will indicate the approximate voltage of the external signal.

The indication on the display may become unstable when no external connection is made and indicates the high sensitivity of the electronics.



For safety reasons the voltage of the external signal must **not exceed** 48 V to earth!

4.3 PC control

If the instrument has been equipped with the optional RS-232 or RS-485 interface, it can be controlled digitally, e.g. from a PC by **PNet peristaltic pump control software**.

This is particularly useful when adding medium to a fermentor or bioreactor to feed the culture which is also growing exponentially (feed pump and harvest pump). Thus, the growth rate and resulting culture activity are considerably increased.



Figure 4-3 Disconnect the pump from the mains. While keeping the direction button ◀▶ pressed, connect the pump to the mains again.



Figure 4-4 The message “A” and two numbers will appear on the display. This number from 00 to 99 is the current address of the pump.



Figure 4-5 To change the address press the buttons ▲ ▲ ▲ under the display until the desired number is obtained.



Figure 4-6 To confirm the address and save it, press the **ON/OFF** button.

5 RECOMMENDATIONS

- ✓ It is advisable to use tubing with a small diameter and high-speed control setting for small flow rates rather than the opposite. This allows a finer selection of flow rates.
- ✓ If possible use the clockwise rotation of the tubing pump. This results in lower friction and the pressure of the liquid is also lower (approx. 0.18 MPa). If a higher pressure is required (approx. 0.2 MPa) use the counter-clockwise rotation.
- ✓ Periodically smear a small amount of petroleum jelly or similar grease on the inner side of the transparent PVC cover of the peristaltic pump. This will increase the reliability and the lifetime of your peristaltic pump. **Do not grease the tubing holding slots.**
- ✓ If, as a consequence of tubing breakage or some other accident, liquid spills into the top of the tubing pump, disconnect the pump from the mains and clean it by removing the liquid and rinsing it with water. The rotor can be completely removed by unscrewing the nut (size M5) on the axis of the rotor and pulling the rotor out by hand or with pliers. After cleaning, grease the axis and replace the rotor by pressing and rotating until the rotor engages on the motor axis.
- ✓ Clean the peristaltic pump with a damp cloth. Mild solvents like ethanol, isopropanol, alkanes are tolerated, if the exposure is short.

If you have any difficulties or questions concerning your MEGAFLOW peristaltic pump, please contact our service office (support@lambda-instruments.com).

6 FOR YOUR SAFETY

Thanks to the use of a plug-in power supply giving only a low voltage of 12 V DC the danger of electrical shock during the use of the MEGAFLOW peristaltic pump has been virtually eliminated, even when an electro conductive solution penetrates the tubing pump.

The peristaltic pump is usually used in a vertical position. The peristaltic pumps can also be stacked, thus allowing optimal use of your expensive laboratory areas.

If the tubing pump is not used for an extended period of time, disconnect it from the mains. A modern miniaturized switching power supply is used, which has only a negligible consumption of electric current when the peristaltic pump is not in use.

7 CONSTRUCTION ADVANTAGES OF THE PUMP

- Instead of small rollers, which are used by most pump producers, LAMBDA peristaltic pumps use ball bearings of a larger diameter with glass beads. This reduces pulsation, friction and mechanical strain on the tubing leading to better pump performance.
- The bearings glide over the tubing so gently, that unlike in other tubing pumps, it is not necessary to prevent the movement of the tubing by any special fixation (stoppers, clamps etc.). As a consequence, the lifetime of the tubing is considerably increased.
- The pressure on the tubing is transmitted gradually through an off-centre lever and spring made of stainless steel. This assures that only minimal pressure is applied to the tubing, which guarantees reliable functioning of the pump without unnecessary deformation of the tubing.
- The spring also reduces the liquid pressure to approx. 2 bar. This is useful when for any reason the line is blocked.
- The asymmetric pump head reduces pulsation and is made from hard, chemically stable material.
- A high quality, Swiss made, BLDC, together with integrated microprocessor electronics assure highest precision of flow rates with no inertia while turning the peristaltic pump on and off.
- The tubing pump dimensions have been minimized. Therefore the MEGAFLOW peristaltic pump is considerably smaller than other products of similar performance. The pump is easy to use and saves expensive laboratory bench space.
- The different remote control options and the possibility of flow integration increase the scope of use of the LAMBDA peristaltic pumps in automatically controlled systems (e.g. in fermentations and cell cultures, chemical synthesis, fraction collection, etc.)

8 UTILIZATION OF LAMBDA PERISTALTIC PUMPS

LAMBDA peristaltic pumps/tubing pumps can be used in the following fields:

- ✓ *Chromatographic techniques:* Liquid chromatography, collecting fractions, sampling, taking samples, gradient elution, pouring gradient gels, gradient formation, etc.
- ✓ *Single use systems:* In contamination sensitive processes, LAMBDA pumps can be used for precise addition or removal of the desired liquid
- ✓ *Fermentation and cell culture techniques:* Addition of nutrients, pH control (buffers, acid, base), antifoam control (adding antifoam agent), feed, harvest, sampling, continuous processes, filling agar plates
- ✓ *Pharmaceutical research:* Process validation, drug trials
- ✓ *Chemical reactions:* Precise addition of liquids for titration, visualization of reaction velocity – hydrolysis of amides, esters, anhydrides, etc.
- ✓ *Biochemical reactions:* Control of oxidation and reduction potential, measurement of enzyme activities, long term reactions



9 TECHNICAL SPECIFICATIONS

9.1 General specification

<i>Type:</i>	LAMBDA MEGAFLOW – microprocessor-controlled programmable peristaltic pump
<i>Programming:</i>	Up to 99 steps of speed and time
<i>Time resolution:</i>	0 to 999 minutes in 1 minute steps 0 to 99.9 minutes in 0.1 minute steps
<i>Accuracy:</i>	± 1%
<i>Reproducibility:</i>	± 0.2 % (electronics)
<i>Tubing:</i>	Silicone tubing or other materials having similar elasticity; inner tubing diameter from 1 to 8 mm and tubing wall thickness of 16pprox.. 2 mm
<i>Flow rate range:</i>	
<i>Minimum:</i>	0.2 ml/hour with 1 mm inner tubing diameter
<i>Maximum:</i>	60,000 ml/hour with 8 mm inner tubing diameter
<i>Non-volatile memory:</i>	Storage of all settings
<i>Maximum pressure:</i>	Approx. 0.18 Mpa in clockwise rotation; approx. 0.2 Mpa in counter-clockwise rotation
<i>Motor:</i>	Microprocessor controlled brushless long life BLDC motor with neodymium magnets
<i>Speed control range:</i>	0 to 999
<i>Interface:</i>	RS-485 (optional) or RS-232 (optional)
<i>Power supply:</i>	90–240 V/50–60 Hz AC plug-in power supply with DC 12V/50W output; possible field operation on 12 V accumulator
<i>Dimensions:</i>	18 (W) x 13 (H) x 16 (D) cm
<i>Weight:</i>	2.5 kg
<i>Safety:</i>	CE, meets IEC 1010/1 norm for laboratory instruments
<i>Operation temperature:</i>	0-40 °C
<i>Operation humidity:</i>	0-90% RH, not condensing
<i>Remote control:</i>	0-10 V; (option 0-20 or 4-20 mA)



For safety reasons the voltage of the external signal must **not exceed** 48 V to earth!

9.2 Remote control (Inputs/outputs)

No.	Colour	Description
1	yellow	(+) input remote speed control 0-10V *)
2	grey	step signal from stepping motor (0 and 12V)
3	green	earth, 0 V
4	brown	+ 12 V
5	white	(+) input remote ON/OFF; 0V = ON, 3–12 V = OFF (this logic can be inversed on demand)
6	pink	earth, ground (GND)
7	red	RS 485 B (-)
8	blue	RS 485 A (+) *) (zero line connected to the contact no. 3)

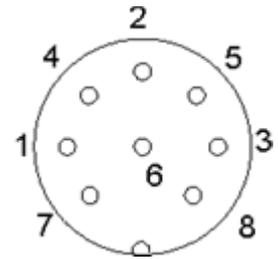


Figure 9-1 8-pole connector

9.3 Input (12 V DC)

Contact No.	Description
1	+ 12 V DC
2	0 V
3	not connected

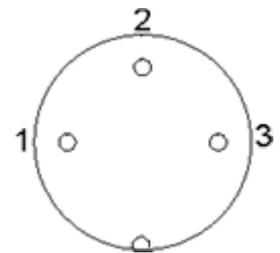


Figure 9-2: 3-pole connector

10 GUARANTEE

LAMBDA provides a two-year guarantee on material and manufacturing defects, if the instrument was used according to the operation manual.

Conditions of guarantee:

- The instrument must be returned with a complete description of the defect or problem. In order to send back the equipment for repair, you will need a returns authorization number from LAMBDA.
- The customer will send the instrument to our service office.
- Damage or loss of items during transport will not be compensated for by LAMBDA.
- Failure to fulfil these requirements will disqualify the customer from compensation.

Serial Number: _____

Guarantee from: _____

11 ACCESSORIES AND SPARE PARTS

11.1 Pump flow INTEGRATOR (Art. No. 4803)

The MEGAFLOW peristaltic pumps and the other LAMBDA instruments like powder dosing instrument LAMBDA DOSER are the only pumps on the market, which allow **a simple and precise integration of the amount of liquid, solid or gas that has been delivered by the pump.**

The electrical impulses, which move the pump motor, are registered and transformed into a direct voltage. This voltage can be measured or recorded by common recorders or voltmeters. The RS 485 interface allows the control

In processes where the pump is controlled, e.g. by a pH-stat during a fermentation or cell culture to keep the pH of the medium constant, it is often important to know when and how much acid or base were added. **This data yield important information about the process, its kinetics and time of completion, etc.**

Another use of the INTEGRATOR is for the **measurement of enzyme activities** (e.g. amidases, esterases, lactamases, lipases, proteases and other enzymes).

The pump-flow INTEGRATOR can now be electronically implemented inside the MULTIFLOW peristaltic and tubing pump and therefore, does not require any additional valuable laboratory bench space.

The activated INTEGRATOR within the LAMBDA pumps allows **new and unusual pump applications** (e.g. gradient making, counter flow elution, liquid chromatography, electronic burette, etc.).

11.2 Pnet control software for peristaltic and syringe pumps, DOSER or MASSFLOW (Art. No. 6600)

Pnet is a PC control software for the remote control of LAMBDA laboratory instruments (peristaltic pumps PRECIFLOW, MULTIFLOW, HIFLOW, MAXIFLOW, MEGAFLOW, syringe pump VIT-FIT, powder dosing instrument DOSER and gas flow controller MASSFLOW).

The pumps are connected to the computer through a RS-232 or RS-485 interface. Up to 6 LAMBDA laboratory instruments and 12 INTEGRATORS can be connected and controlled simultaneously.

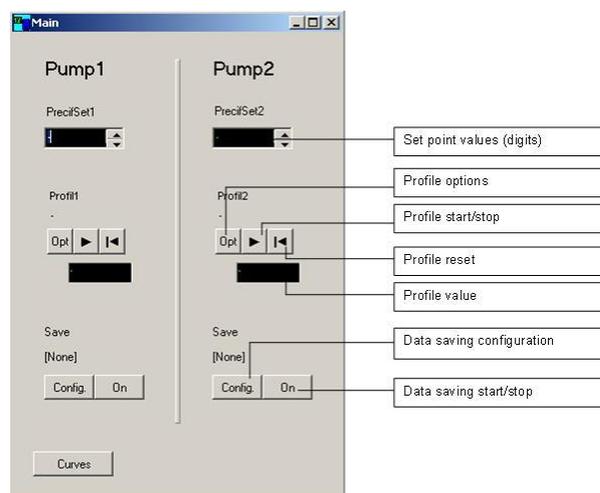


Figure 11-1 Remote control of the peristaltic pump with the PNet control software

11.3 List of accessories and spare parts

Art. No.	Accessories
4803	<i>PUMP-FLOW INTEGRATOR (for LAMBDA pumps, DOSER and MASSFLOW)</i>
4810	<i>Pump remote control (analog and digital) cable, 8 poles (open ends)</i>
4802	<i>Pump ON/OFF remote control cable, 2 poles (open ends)</i>
4823	<i>Foot switch for ON/OFF control, 8 poles</i>
4823-b	<i>Foot switch for ON/OFF control (bistable switch)</i>
4824	<i>Cable for inverted analog ON/OFF control, 8 poles</i>
 Interface and control software	
4822	<i>RS232 interface (for connection of the instruments to the serial port)</i>
4816	<i>RS485 interface (for connection of the instruments to the serial port)</i>
4817	<i>RS232/485 converter</i>
4817-kit	<i>RS485 connection kit (for connection to a serial port or USB port)</i>
4818	<i>Power supply for RS232/485 converter (5V/1W)</i>
4819	<i>RS-line connection cable (serial)</i>
6600	<i>PNet control software for peristaltic and syringe pumps, DOSER or MASSFLOW</i>
800202	<i>Quadruple plug box (Power and RS-connection for up to 4 LAMBDA laboratory instruments)</i>
 Spare parts	
6003	<i>Plug-in power supply (12V/65W) for MEGAFLOW (Plug type: AU, EU, UK, US)</i>
4805	<i>Roller</i>
4806	<i>Stainless steel spring</i>
4807	<i>Eccentric lever</i>
4808	<i>Rotor</i>
4809	<i>Cover</i>
4811	<i>Pump head</i>
800100-26	<i>Silicone tubing 2/6 mm, 5 m</i>
800100-36	<i>Silicone tubing 3/6 mm, 5 m</i>
800100-48	<i>Silicone tubing 4/8 mm, 5 m</i>
800100-610	<i>Silicone tubing 6/10 mm, 5 m</i>
800100-812	<i>Silicone tubing 8/12 mm, 5 m</i>
800113	<i>Stainless steel tubing clamp</i>

12 APPENDIX

12.1 RS communication protocol for LAMBDA PERISTALTIC and SYRINGE PUMPS, powder DOSER and MASSFLOW gas flow regulator

12.1.1 Format of data sent by the PC to the pump and back

Data sent by the PC: #ss mm a ddd qs c

Data sent back by the pump: <mm ss a ddd qs c

where,

- #** is the first sign of a command sent by PC
- <** is the first sign of a message sent by pump
- ss** is the address of the pump
- mm** is the address of the PC
- a** is the command for the sense of rotation:
- r** for clockwise (cw) rotation (to the right)
- l** for counter-clockwise (ccw) rotation (to the left) (not for DOSER and MASSFLOW)
- ddd** is the speed of rotation (3 ASCII numbers from 0 to 9; sent from the highest order digit to the lowest order digit)
- qs** is the control sum in HEX format (2 ASCII signs of the type 0...9ABCDEF)
- c** is the end sign cr (carriage return) The pump will fulfill the task and block any manual command on the pump front panel.

12.1.2 Commands not containing data

- # ss mm **g** qs c activates the local command of the pump
- # ss mm **s** qs c the pump is stopped
- # ss mm **G** qs c to send pump data to the PC

12.1.3 Checksum control

The PC sends: #0201r123EEcr

The control sum (checksum) qs is made in the following way (only the **last byte** (2 ASCII characters of the type 0...9ABCDEF) is taken):

#	0	2	0	1	r	1	2	3	EE (last byte)	cr
23h	+30h	+32h	+30h	+31h	+72h	+31h	+32h	+33h	=1EEh	0Dh

12.1.4 Format of the data transmission

Speed: 2400 Bd (Baud)
8 data bits, odd parity, 1 stop bit

12.2 Examples

Address of the PC: 01
Address of the pump: 02

The PC sends: #0201r123EEcr
The pump will rotate cw at the speed of 123

The PC sends: #0201G2Dcr
The answer of the pump: <0102r12307cr

The PC sends: #0201l123E8cr
The pump will rotate ccw at the speed of 123. (not for DOSER and MASSFLOW)

The PC sends: #0201s59cr
The pump stops.

The PC sends: #0201g4Dcr
The pump will go to the local command (pump front panel is activated).

12.3 How to set the MEGAFLOW peristaltic pump address?

To look up/modify the instrument address, disconnect the MEGAFLOW peristaltic pump from mains. (Refer figures 4-3, 4-4, 4-5, 4-6 under the section 4.3 PC control)

Press the ◀▶ button continuously and at the same time connect the MAXIFLOW peristaltic pump to the mains again. The message “A” and two numbers will appear on the display. This number from 00 to 99 is the current address of the peristaltic pump.

To change the address, press the buttons ▲▲▲ under the display until the desired number is obtained.

To confirm and save the address, press the ON/OFF button.

12.4 RS-connection scheme

The 8-pole DIN connector “REMOTE” is used for the remote control and RS-485 connection. When the optional RS-485 interface is available the pins are used as follows:

No.	Colour	Description
1	yellow	(+) input remote speed control 0-10V *)
2	grey	step signal from stepping motor (0 and 12V)
3	green	earth, 0 V
4	brown	+ 12 V
5	white	(+) input remote ON/OFF; 0V = ON, 3–12 V = OFF (this logic can be inversed on demand)
6	pink	earth, ground (GND)
7	red	RS 485 B (-)
8	blue	RS 485 A (+)

*) (zero line connected to the contact no. 3)

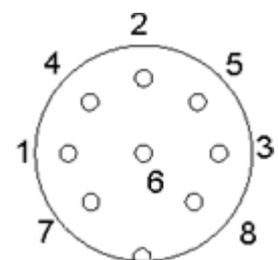


Figure 12-1: 8 pole connector

12.5 RS communication protocol for the on-board INTEGRATOR (optional)

12.5.1 Communication between the PC and the INTEGRATOR of the LAMBDA instrument

From the PC to the INTEGRATOR:

#ss mm z qs c

From the INTEGRATOR to the PC:

<mm ss = qs c confirmation of the reception of a command
<mm ss dddd qs c sending of the requested data

where,

- #** is the first sign of a command sent by the MASTER (PC)
- <** is the first sign of a message sent by the SLAVE (INTEGRATOR)
- ss** is the address of the subordinate station (address of the instrument with integrated INTEGRATOR)
- mm** is the address of the commanding station (PC)
- z** is a command (see below): small letters indicate a command, capital letters request data transfer from the subordinate station
- =** confirmation of reception
- aa** new address of the subordinate station (ss) (two numbers and possibly other ASCII characters A B C D E F)
- dddd** transferred data (values are two bytes in hexadecimal form. Single bytes are transformed into two ASCII characters 0,...,9,A,B,C,D,E,F)
- qs** is the control sum (obtained by the addition modulo 256 of binary values of all preceding characters including the leading sign) in HEX format (2 ASCII signs of the type 0...9ABCDEF)
- c** is the end sign cr (carriage return)

12.5.2 Commands for the INTEGRATOR

- n** reset (sets the INTEGRATOR to zero)
- i** start of integration
- e** stop of integration
- I** sends the integrated value
- N** sends the integrated value and sets the integrator to zero
- L** sends the integrated value ccw rotation (to the left) (not for DOSER)
- R** sends the integrated value of cw rotation (to the right)

12.5.3 Examples

Address of the PC: 01
 Address of the instrument with on-board INTEGRATOR: 02

The PC sends: #0201I2Fcr

The control sum (checksum) qs is made in the following way (only the **last byte** (2 ASCII characters of the type 0...9ABCDEF) is taken):

#	0	2	0	1	I	2F (last byte)	cr
23h	+30h	+32h	+30h	+31h	+49h	=12Fh	0Dh

The PC sends: #0201i4Fcr
 i.e. in hexadecimal form: 23h 30h 32h 30h 31h 69h 34h 46h 0Dh
 This means: For a subordinate station (SLAVE) with address 02 from commanding station (MASTER) with address 01
 Start of integration
 The control sum is 14Fh (last byte: **4F**); end of message cr (carriage return)
 The INTEGRATOR answers: <0102=3Ccr

The PC sends: #0201N34cr
 The INTEGRATOR answers: <0102N03C225cr (integrated value is 03C2h)
 and resets to zero

The PC sends: #0201e4Bcr
 The integration will be stopped and the command will be confirmed.
 The INTEGRATOR answers: <0102=3Ccr



LAMBDA Laboratory Instruments

Sihlbruggstrasse 105
 CH-6340 Baar
 SWITZERLAND – EUROPE
 Tel.: +41 444 50 20 71
 Fax: +41 444 50 20 72

E-mail: support@lambda-instruments.com
 Web: www.lambda-instruments.com

LAMBDA CZ s.r.o.

Lozibky 1
 CZ-61400 Brno
 CZECH REPUBLIC – EUROPE

Hotline: +420 603 274 677

www.peristalticpump.info
www.peristaltic-pumps.eu